

Introduction

$\gamma\delta$ T cells are promising effector cells for adoptive immunotherapy. Thanks to their MHC-independent antigen recognition and the detection of stress ligands^{1,2}, they are suitable for allogeneic clinical interventions. V γ 9V δ 2 T cells, the main population in human peripheral blood, have the advantage of naturally recognizing malignant cells, especially those of hematopoietic malignancies³ by different mechanisms, such as sensing the upregulation of the mevalonate pathway happening in cancer cells⁴. This mechanism can be used to successfully expand these cells with


Zoledronate and cytokines⁵. Moreover, their safety has already been demonstrated in both autologous^{6,7} and allogeneic⁸ clinical trials for hematological and solid tumors. However, genetic modification of $\gamma\delta$ T cells has been shown to be challenging due to inefficient transduction. To enhance their natural cytotoxicity, a suitable engineering method is thus needed. Here, we have established an easily scalable protocol to produce CAR V γ 9V δ 2 T cells with a high transduction efficiency combined with a stable final phenotype.

Results

1 $\gamma\delta$ T cells are efficiently expanded with a feeder-cell free method


We activated and expanded $\gamma\delta$ T cells from peripheral blood mononuclear cells (PBMC) using Zoledronic acid 5 μ M, IL-2 100 IU/ml and IL-15 100 IU/ml in TexMACSTM medium in the presence of human AB serum. We monitored the cellular composition, the pheno-

type and the activation profile by flow cytometry during two weeks of culture. We expanded V γ 9V δ 2 T cells from 5 donors at least 1000 times, reaching a purity above 75%. The final product consisted largely of $\gamma\delta$ T cells and NK cells (Figure 1).

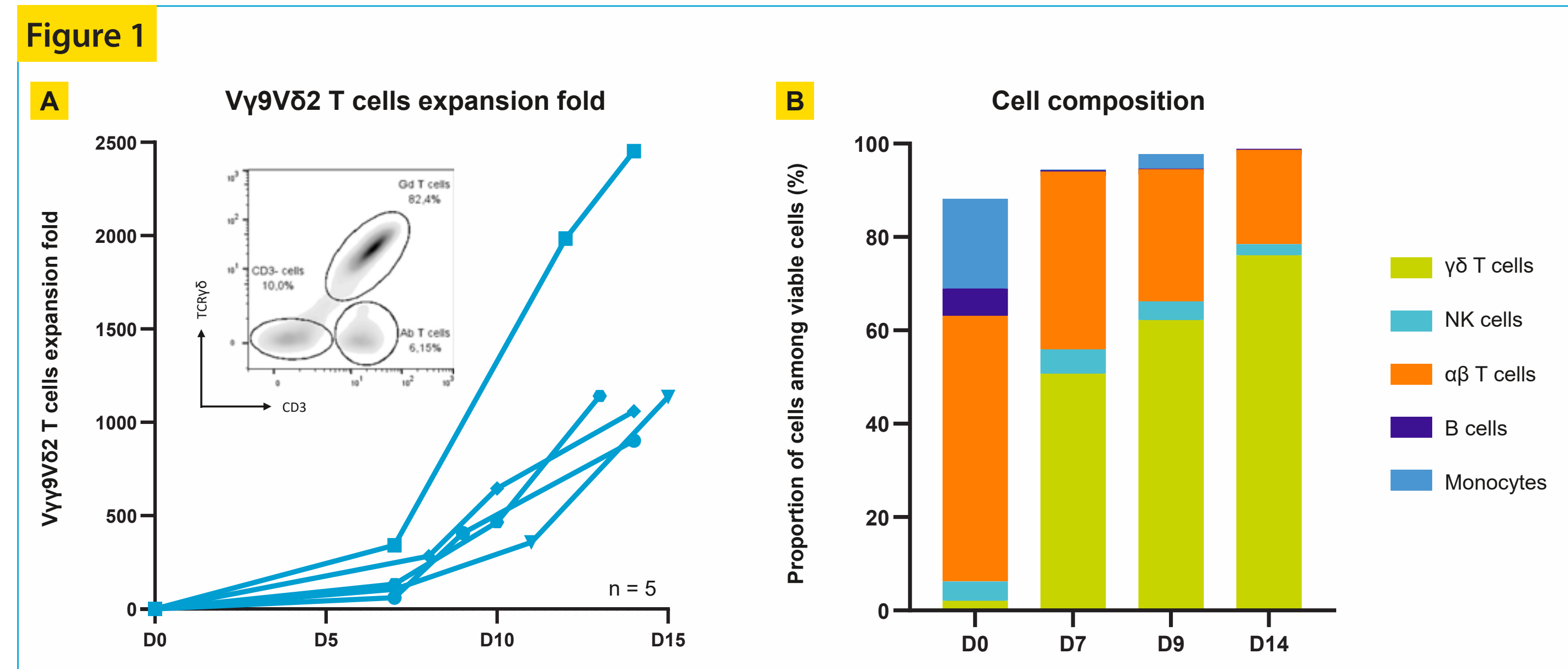


In all donors, the $\gamma\delta$ T cells showed a memory phenotype (Figure 2A), and a high expression of the ac-

tivation markers CD69, CD56 and HLA-DR at the end of the cultivation process (Figure 2B).

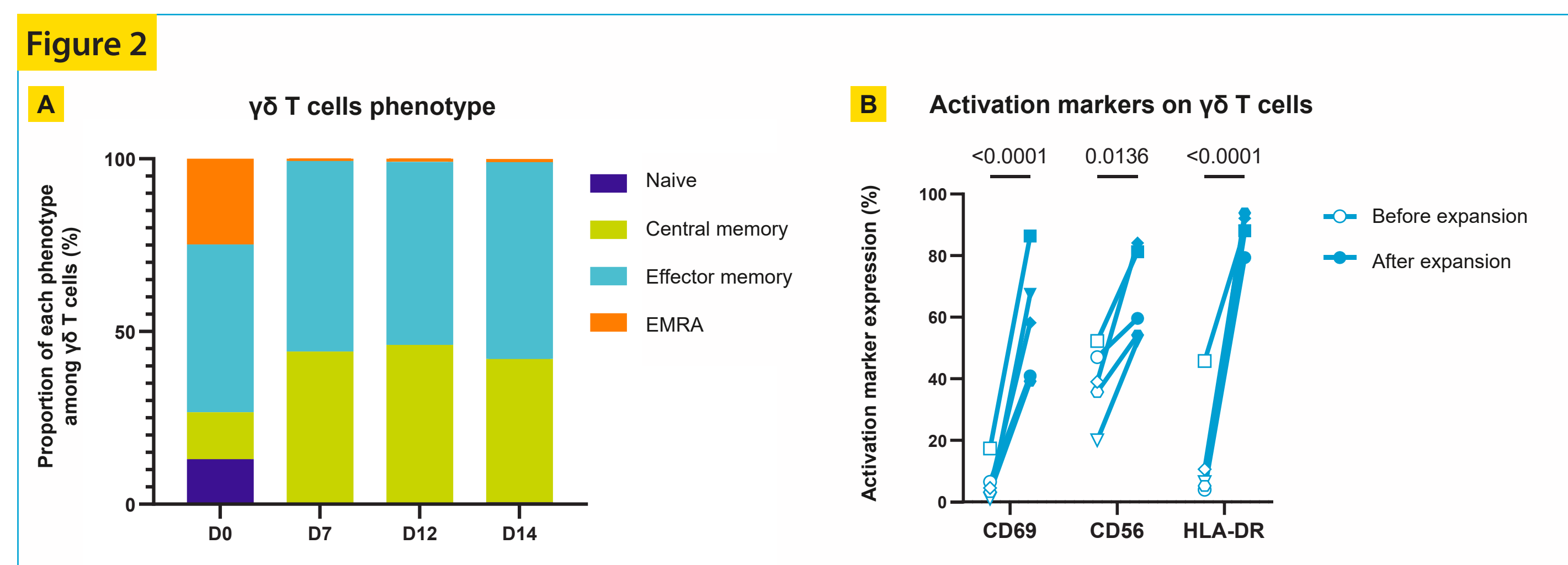

We performed $\alpha\beta$ T-cell depletion using magnetic beads. The depletion was then performed with LS columns (Figure 3A) or on the CliniMACS Prodigy (Figure 3B).


A depletion of at least 4 logs was reached with both methods. $\gamma\delta$ T cells could then be expanded from the PBMCs.

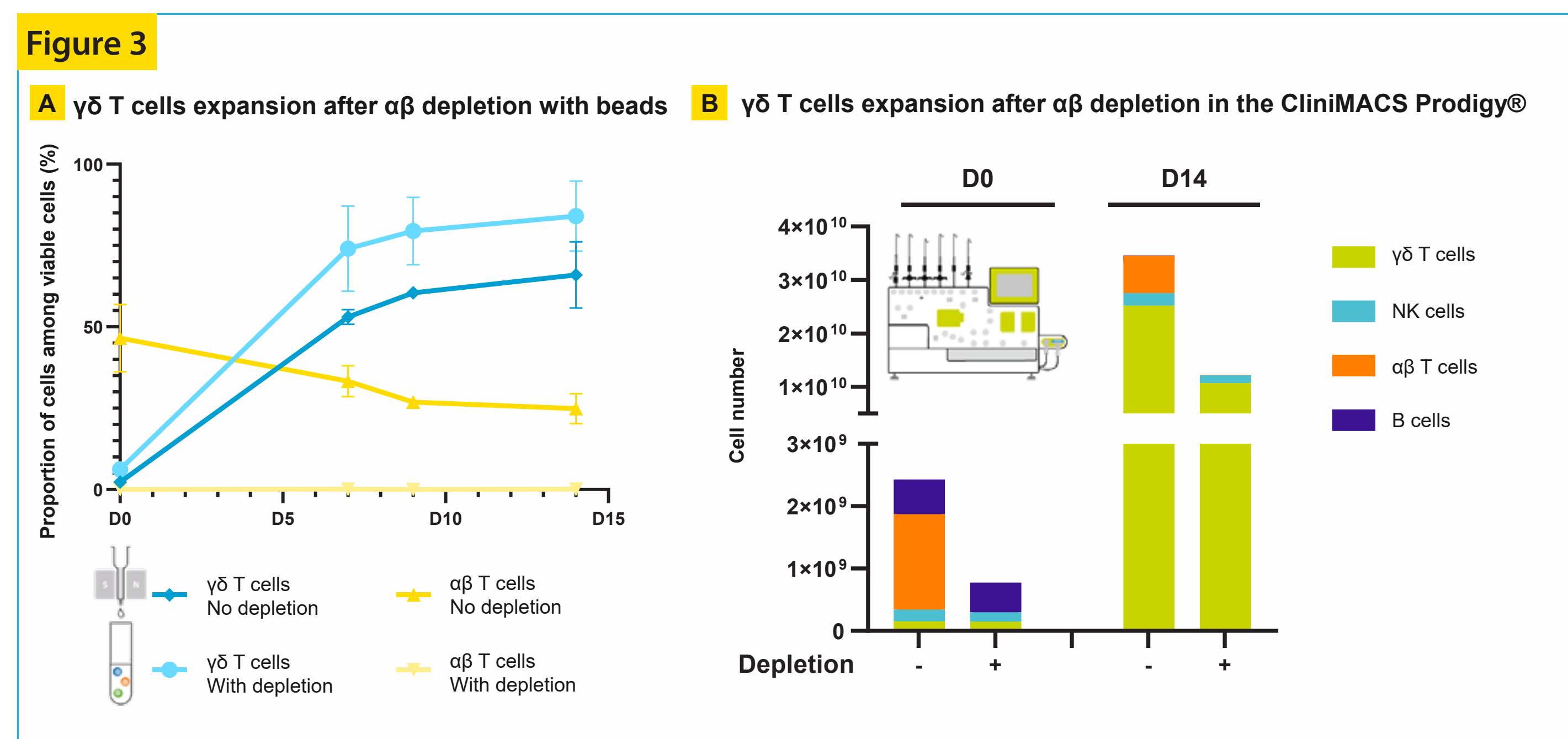

2 BaEV-pseudotyped lentiviral vector yields a high transduction efficiency of $\gamma\delta$ T cells at a low MOI

We compared the VSV-G pseudotype with BaEV at a MOI of 0.5 with CD19 CAR and CD33 CAR constructs. Vectofusin[®]-1 was used at a concentration of 10 μ g/ml to enhance the transduction efficiency of the BaEV pseudotype. Vectofusin-1[™] did not increase transduction efficiency with VSV-G (data not shown). CAR expression and the number of $\gamma\delta$ T cells were determined 7 days after transduction. The combination of


the BaEV pseudotype and Vectofusin[®]-1 resulted in a 86.1% \pm 9.4% transduction efficiency, which was significantly superior to VSV-G transduction efficiency for both CAR constructs tested (Figure 4A, $p = 0.0006$). The median fluorescence intensity was similarly higher with the BaEV pseudotype. The transduction did not impact the $\gamma\delta$ T cell expansion (Figure 4B).

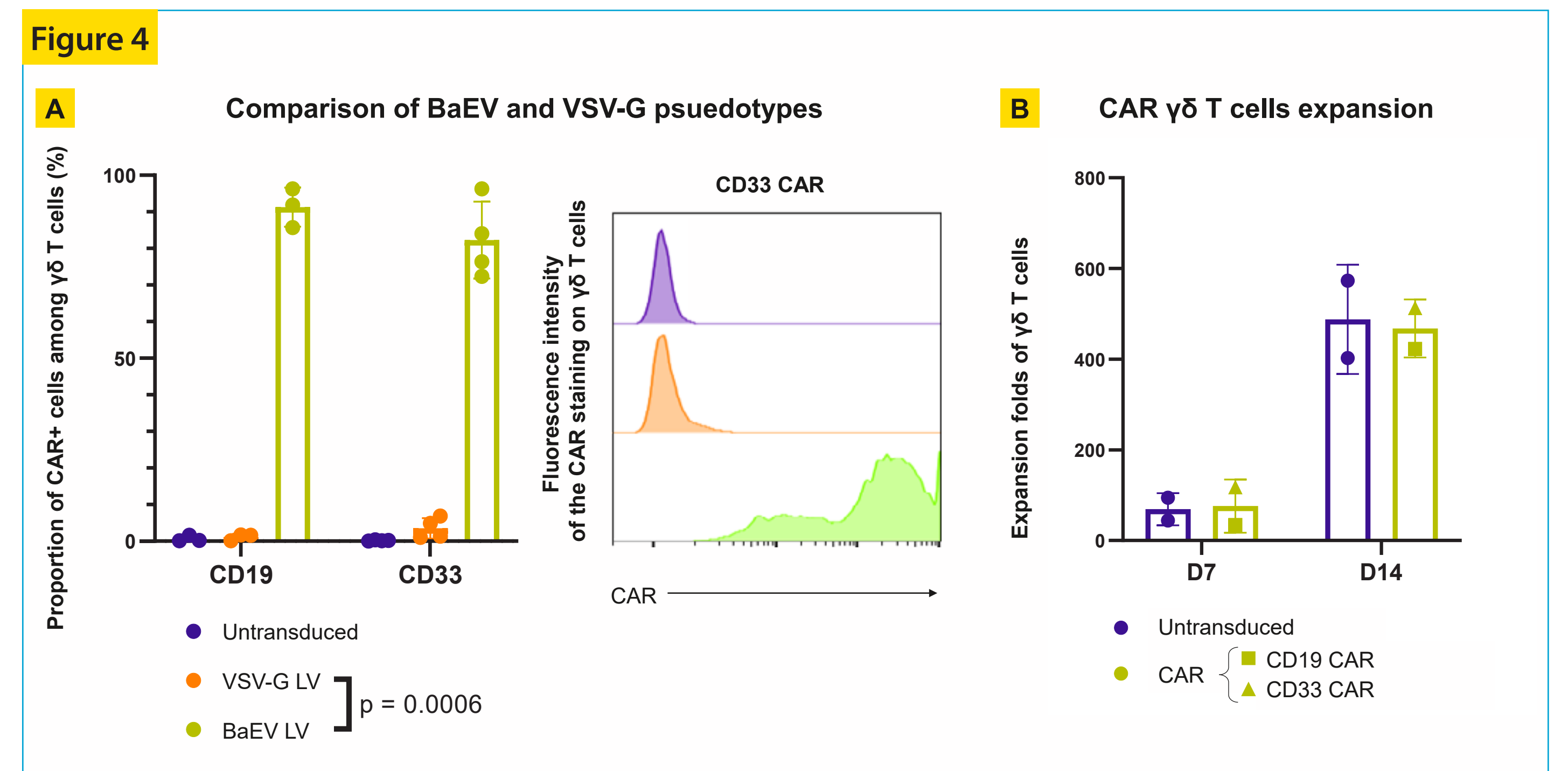

3 CAR $\gamma\delta$ T cells efficiently kill hematological tumor cell lines *in vitro*

We tested the killing efficiency of $\gamma\delta$ T cells after $\alpha\beta$ T-cell depletion, transduction and further expansion with our protocol. The untransduced and CD19 or CD33 CAR $\gamma\delta$ T cells were co-cultured with RS4-11 or OCI-AML3 WT and OCI-AML3 CD33 KO tumoral cell lines expressing GFP and Luciferase reporters. After 4h of co-culture, the Luciferase expression was measured with the ONE-Glo[™] Luciferase Assay System (Promega) as a marker of the number of target cells.

The growth of the target cells was monitored in parallel in the Incucyte[®] S3 System (Sartorius). CD19 CAR $\gamma\delta$ T cells, but not untransduced $\gamma\delta$ T cells, efficiently killed RS4-11 cells in both systems (Figure 5A). Similarly, CD33 CAR $\gamma\delta$ T cells killed OCI-AML3 WT cells better than untransduced $\gamma\delta$ T cells. The superior killing was only observed with the CD33 CAR and not with CD33 KO cells (Figure 5B).

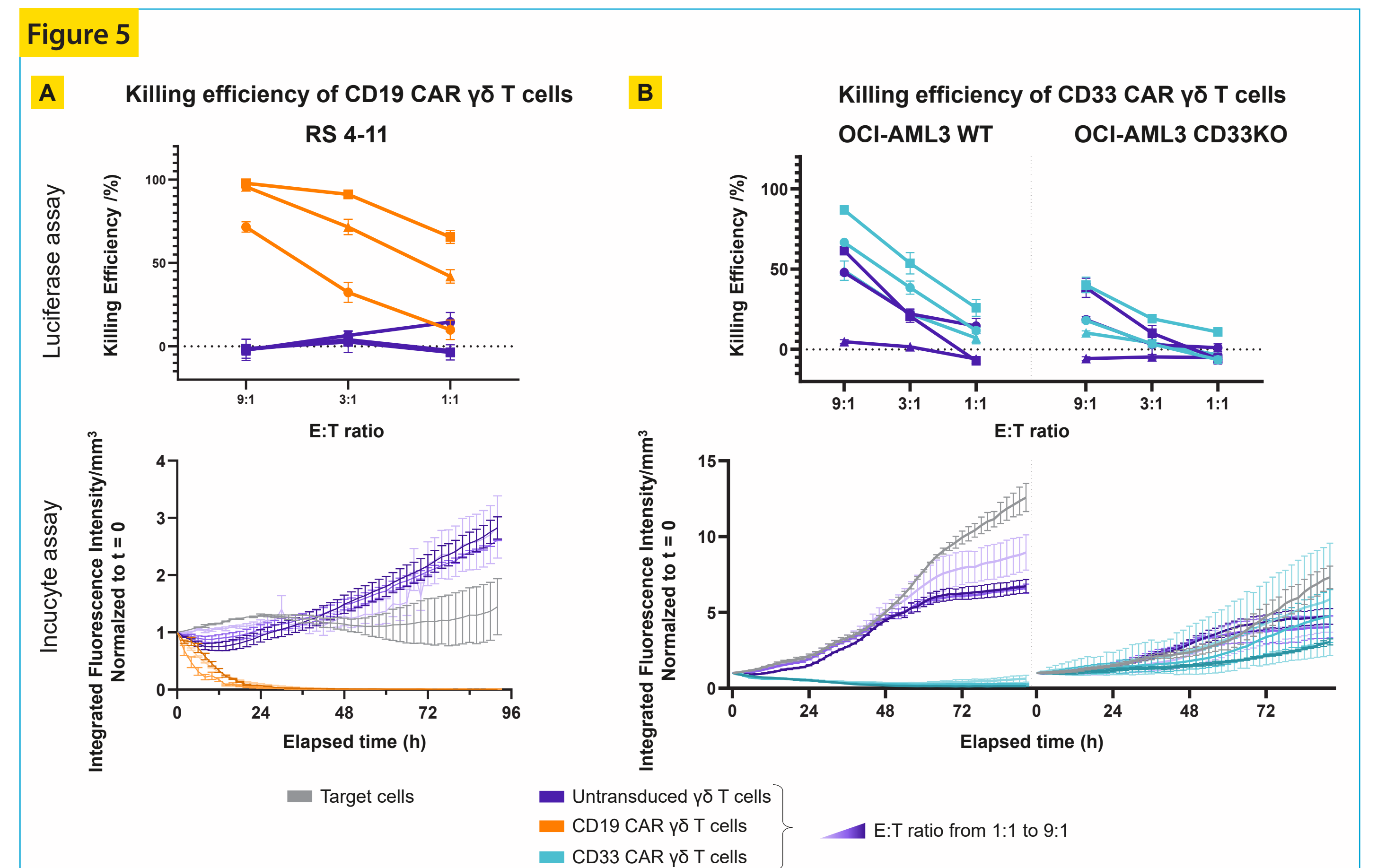

Conclusion

- V γ 9V δ 2 T cells can successfully be expanded by a feeder cell-free protocol after $\alpha\beta$ T cells depletion.
- The BaEV pseudotype combined with Vectofusin[®]-1 results in a high CAR transduction efficiency, without loss of expansion potential.
- Transducing V γ 9V δ 2 with CAR constructs enhances their natural cytotoxicity.
- This method is suitable for clinical production in a closed system such as the CliniMACS Prodigy[®].

References

1. Groh, V. et al. Recognition of stress-induced MHC molecules by intestinal epithelial gammadelta T cells. *Science* 279, 1737–1740 (1998).
2. Hayday, A. C. $\gamma\delta$ cells: a right time and a right place for a conserved third way of protection. *Annu. Rev. Immunol.* 18, 975–1026 (2000).
3. Fisch, P. et al. Control of B cell lymphoma recognition via natural killer inhibitory receptors implies a role for human V γ 9V δ 2 T cells in tumor immunity. *Eur. J. Immunol.* 27, 3368–3379 (1997).
4. Constant, P. et al. Stimulation of Human $\gamma\delta$ T Cells by Nonpeptide Mycobacterial Ligands. *Science* 264, 267–270 (1994).
5. Künzmann, V. et al. Stimulation of $\gamma\delta$ T cells by aminobisphosphonates and induction of antiplasma cell activity in multiple myeloma. *Blood* 96, 384–392 (2000).
6. Benounna, J. et al. Phase-I study of Imnacel gammadelta, an autologous cell-therapy product highly enriched in gamma9delta2 T lymphocytes, in combination with IL-2, in patients with metastatic renal cell carcinoma. *Cancer Immunol. Immunother.* 57, 1599–1609 (2008).
7. Kakimi, K. et al. Adoptive transfer of zoledronate-expanded autologous V γ 9V δ 2 T-cells in patients with treatment-refractory non-small-cell lung cancer: a multicenter, open-label, single-arm, phase 2 study. *J. Immunother. Cancer* 8, e001185 (2020).
8. Xu, Y. et al. Allogeneic V γ 9V δ 2 T-cell immunotherapy exhibits promising clinical safety and prolongs the survival of patients with late-stage lung or liver cancer. *Cell. Mol. Immunol.* (2020).